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GymBeam

SAFETY WARNING

Before using this product, we recommend that you carefully study the enclosed manual and keep it in case of need in future.

1. The product is intended for external use only (on the skin only).
2. The product is not a substitute for an orthopaedic brace and is not intended for therapeutic purposes.
3. Do not overheat in the microwave as this may cause the gel pack to burst or leak.
4. Do not use the product for more than 20 minutes. Reapply every 1–1.5 hours if needed.
5. The gel pack shouldn't come into direct contact with your skin as it may cause hot or cold burns. Always use a brace as a protective layer.
6. Do not puncture or tear the gel pack. Protect it from contact with sharp objects.
7. Check the gel pack for punctures before each use. If there's a puncture, throw the gel pack away and wash your hands with soap.
8. The gel pack is non-toxic and latex free. If the gel pack substance comes to contact with your skin or eyes, immediately rinse the affected area with warm water. If the gel pack substance is swallowed, please consult a physician immediately.
9. In case of swelling, discomfort or pain, stop using the product immediately. If pain or discomfort persists, consult your physician.
10. Adult supervision is recommended when using this product on children.

SPECIFICATION

Therapeutic Ankle & Foot Brace Hot-Cold is designed to be used for recovery period after injuries or surgery, or for any other ankle/foot problems. Utilizes hot/cold therapy packs.

PRODUCT DETAILS

Big gel pack size: 32.7x13.8 cm; **Small gel pack size:** 16.5x10.4 cm; **Big gel pack weight:** 220 g; **Small gel pack weight:** 130 g; **Gel pack material:** PVC & NYLON; **Brace material**: Neoprene; **Brace size:** 45.7x29.5 cm; **Brace weight:** 112 g

CARE INSTRUCTIONS

Remove gel packs before washing. Wash after every 3–4 uses in cold water with antibacterial soap, rinse thoroughly, squeeze excess moisture and air dry. Do not machine wash. Do not bleach.

COOLING INSTRUCTIONS

To cool gel packs remove them from the brace. Place them flat in the freezer for at least 2 hours. Do not use the product for more than 20 minutes. Make sure that the gel pack doesn't come into direct contact with your skin as it may cause cold burns. Reapply every 1–1.5 hours if needed.

HEATING INSTRUCTIONS

Remove gel packs from the brace and place them in the microwave. Heat for 20 seconds (600-800W) and if gel packs are not warm enough, microwave for 10 more seconds to get the desired temperature. Do not overheat and do not wear if uncomfortably hot. Make sure that the gel pack doesn't come into direct contact with your skin as it may cause hot burns. If the gel pack is too hot, let it cool before application. Do not use the product for more than 20 minutes. Reapply every 1–1.5 hours if needed.


HOW TO USE THE ANKLE & FOOT BRACE


Step 1

Insert heated or cold gel packs into the pockets on the inside of the brace as illustrated on the picture. There are 2 pockets.

Step 2

Place the brace over the ankle. Use Velcro straps to secure the brace in place.

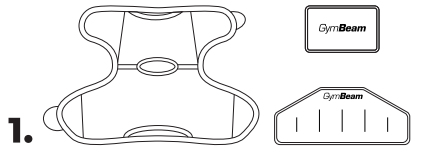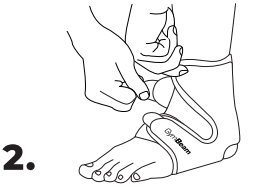WARRANTY

The warranty is valid for 2 years from the date of purchase and will be void in the following cases:

- Damage caused by accident, misuse, improper handling or transportation.
- Unauthorized repairs or tampering with the product.
- Failure to use the product according to the provided instructions.
- Using the product for purposes other than its intended use.
- Failure to follow maintenance, operation, and service instructions.
- Damage exceeding the value of the product.
- Natural wear and tear of product components resulting from normal use.
- Exposure to adverse external factors such as sunlight or other radiation, liquid ingress, or electrostatic discharge.
- Damage to the product due to natural disasters, unauthorized intervention, or mechanical damage caused by the buyer (e.g., during transportation, using improper cleaning agents).
- Deterioration of the product's condition due to improper storage, handling, or care by the user.

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